KERN BALANCES & TEST SERVICES CATALOGUE 2021

Hanging scale KERN HDB-N



Practical hanging scale for rapid, mobile weighing

Features

- With the TÜV certification mark, the scales meet the requirements of the standard EN 13155 (Non-fixed load lifting attachments/ breakage resistance) and EN 61010-1 (Electrical safety)
- Ideal for rapid control in goods-in and goods-out
- Also essential in the private sector to determine the weight of fish, game, fruits, bicycle parts, suitcases etc.
- Hold function: For easy reading of the weighing result, the display can be "frozen" in different ways. Either automatically when the weighing value remains unchanged or manually by pressing the Hold key
- Hook, steel, can be hinged

STANDARD

Peak load display (peak hold)

Technical data

- · LCD display, digit height 12 mm
- Ready for use: Batteries included, 2×1.5 V AAA, operating time up to 180 h
- Further weighing units: kg, lb, N
- Permissible ambient temperature 5 °C/35 °C

Accessories

• Tare pan with pan beam, details see Accessories, KERN CH-A01N

CAL EXT UNIT MOVE	BATT DMS 1 DAY +3	DAYS		
Model	Weighing capacity	Readability	Dimensions ho W×D×H	
	[] [] [] [] [] [] [] [] [] []	1 .1	W^D^H	

OPTION

DAkkS

Model	Weighing capacity	Readability	Dimensions housing	Net weight	Option	
			W×D×H		DAkkS Calibr. Certificate	
	[Max]	[d]		approx.	DAkkS	
KERN	kg	g	mm	kg	KERN	
HDB 5K5N	5	5	72×25×105	0,20	963-127H	
HDB 10K10N	10	10	72×25×105	0,20	963-128H	

S TUV SUD Production monitored Type tested





KERN BALANCES & TEST SERVICES CATALOGUE 2021

KCP

PROTOCOL

GLP

INTERN

PRINTER

PCS

RECIPE

RECIPE

- 88'

SUM

PERCENT

C

UNIT

- → +<

TOL

^-

digital systems GLP/ISO log:

connection GLP/ISO log:

printers

Piece counting:

Recipe level A:

Recipe level B:

Totalising level A:

value (100 %)

Weighing units:

Hold function:

KERN Communication Protocol (KCP):

It is a standardized interface command set for

KERN balances and other instruments, which

devices featuring KCP are thus easily integrated

with computers, industrial controllers and other

The balance displays serial number, user ID,

With weight, date and time. Only with KERN

Reference quantities selectable. Display can

The weights of the recipe ingredients can

be added together and the total weight of

Internal memory for complete recipes with

The weights of similar items can be added

Determining the deviation in % from the target

Can be switched to e.g. nonmetric units at the

(Checkweighing) Upper and lower limiting can

be programmed individually, e.g. for sorting and

dosing. The process is supported by an audible

(Animal weighing program) When the weighing

conditions are unstable, a stable weight is calculated as an average value

or visual signal, see the relevant model

touch of a key. See balance model. Please refer

together and the total can be printed out

name and target value of the recipe ingredients.

be switched from piece to weight

the recipe can be printed out

User guidance through display

Percentage determination:

to KERN's website for more details

Weighing with tolerance range:

weight, date and time, regardless of a printer

allows retrieving and controlling all relevant parameters and functions of the device. KERN



Pictograms



Internal adjusting: Quick setting up of the balance's accuracy with



Adjusting program CAL:

For quick setting up of the balance's accuracy. External adjusting weight required

internal adjusting weight (motordriven)



Easy Touch: Suitable for the connection, data transmission and control through PC, tablet or smartphone.



Memory:

Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



Alibi memory:

Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard

Data interface RS-232:

• 6550.• To connect the balance to a printer, PC or RS 232 network



RS-485 data interface:

To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible



USB data interface:

To connect the balance to a printer, PC or other peripherals

Bluetooth* data interface:

To transfer data from the balance to a printer, PC or other peripherals



*

WiFi data interface:

To transfer data from the balance to a printer, PC or other peripherals



Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



Analogue interface:

to connect a suitable peripheral device for analogue processing of the measurements



Interface for second balance:

For direct connection of a second balance



Network interface:

For connecting the scale to an Ethernet network






water splashes IPxx: The type of protection is shown in the pictogram

Protection against dust and

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners

UNDER the balance

Ę.

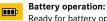



Ready for battery operation. The battery type

Suspended weighing:



is specified for each device

Load support with hook on the underside of



Rechargeable battery pack: Rechargeable set

Universal mains adapter:

with universal input and optional input socket MULTI adapters for A) EU, CH, GB; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS



Mains adapter:

230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available

Power supply:



Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request



Weighing principle: Strain gauges:

Electrical resistor on an elastic deforming body



Weighing principle: Tuning fork:

A resonating body is electromagnetically excited, causing it to oscillate



Weighing principle: Electromagnetic force compensation:

Coil inside a permanent magnet. For the most accurate weighings



Weighing principle: Single cell technology:



Advanced version of the force compensation principle with the highest level of precision



The time required for verification is specified +3 DAYS in the pictogram

DAkkS calibration possible (DKD): DAkkS The time required for DAkkS calibration is +3 DAYS shown in days in the pictogram

Factory calibration (ISO):



The time required for Factory calibration is shown in days in the pictogram



Package shipment:



The time required for internal shipping preparations is shown in days in the pictogram

Pallet shipment:



Your KERN specialist dealer:

The time required for internal shipping preparations is shown in days in the pictogram

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAkkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAkkS calibration laboratory today is one of the most modern and bestequipped DAkkS calibration laboratories for balances, test weights and force-measurement in Europe

Thanks to the high level of automation, we can carry out DAkkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- · DAkkS calibration of balances with a maximum load of up to 50 t
- · DAkkS calibration of weights in the range of 1 mg 2500 kg · Volume determination and measuring of magnetic susceptibility (magnetic
- characteristics) for test weights · Database supported management of checking equipment and reminder service
- · Calibration of force-measuring devices
- · DAkkS calibration certificates in the following languages DE, EN, FR, IT, ES, NL, PL
- · Conformity evaluation and reverification of balances and test weights